

DORAI

SAFETY & RELIABILITY

Your endoscopes dried,
safe & always ready to use

MEDICAL

HEALTH
CARE

NANO

SCIENCE

SERVICE



cbcmind.com



HIGH CAPACITY & FLEXIBILITY

Vertical storage capacity up to 12 endoscopes in a small volume.


Internal sliding walls offer quick and easy access to the endoscopes.

Available in two versions: with single door or pass through.



COMPLETE DRYING

Filtered air through a HEPA filter H14 and twelve independent micro-compressors allow the drying of both the external surface and each inner endoscope's channel.

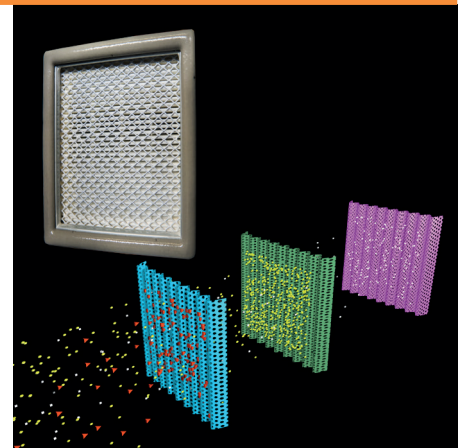


ACCESSIBILITY

A glass door and LED lamps give a quick overview of all available endoscopes.



MAXIMUM SAFETY

Various sensors control all cabinet's parameters and endoscopes' conditions reporting any potential faults.

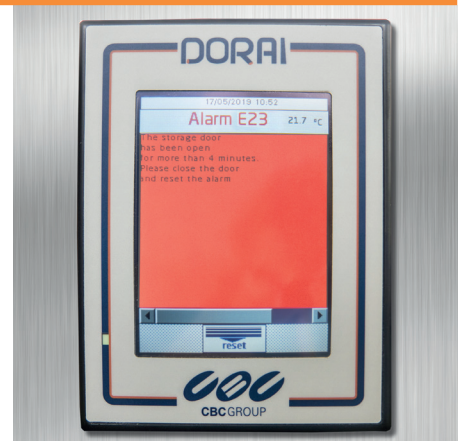

ENDOSCOPE CONNECTION SETS

Specific endoscope connection sets guarantee a perfect drying of the inner channels.



USER FRIENDLY INTERFACE

An intuitive touch-screen display helps to identify the endoscopes' status.



COMPLETE TRACEABILITY

Barcode scanner to easily identify operators and scopes.

Reports can be saved and printed offering a safe and complete traceability.



STANDARDS COMPLIANCE

Compliant with EN 16442:2015 "Controlled environment storage cabinet for disinfected thermolabile endoscopes".

All the tests have been performed by independent laboratory.





SPECIFICATIONS AND TECHNICAL DATA

Weight	kg	350		
Width / Depth / Height *Pass through version	mm	1100 / 860 (920*) / 2100		
Noise level	dB(A)	<50		
Air feed	m ³ /h	100		
Ambient temperature	°C	Min 5°C Max 35°C		
Relative humidity (max.)	HR%	90		
Voltage	V	230	230	120
Frequency	Hz	50	60	60
Power consumption (max.)	W	700	700	700
Current	A	3	3	5.8
Main switch	A	10	10	10

Specifications, design, and accessories are subject to change without any notice or obligation by the manufacturer.